



Ready-to-Wear Medical LegWear - Selection Chart

	Name and Style	S		M		L		L		XL		XL		Compression	Colour						
		Open Toe	Closed Toe	Open Toe	Closed Toe	Open Toe	Closed Toe	Open Toe	Closed Toe	Open Toe	Closed Toe	Open Toe	Closed Toe								
								Full Calf	Full Calf			Full Calf	Full Calf								
FORMAL	UltraSheer Stockings Knee High	75165-30	75165-25	75166-30	75166-25	75167-30	75167-25			75168-30	75168-25			15-20 mmHg	● Natural						
			75165-04		75166-04		75167-04				75168-04				● Black						
			73633-00		73633-01		73633-02				73633-03				● Diamond Pattern Black						
			75201-19		75202-19		75203-19				75204-19				● Natural						
			75201-04		75202-04		75203-04				75204-04			● Black							
			73634-00		73634-01		73634-02				73634-03			● Diamond Pattern Black							
			75241-21		75242-21		75243-21				75244-21			● Natural							
			75241-04		75242-04		75243-04				75244-04			● Black							
	UltraSheer Stockings Thigh High		75169-22		75170-22		75171-22				75172-22			15-20 mmHg	● Natural						
			75169-04		75170-04		75171-04				75172-04				● Black						
			73635-00		73635-01		73635-02				73635-03				● Diamond Pattern Black						
			75205-23		75206-23		75207-23				75208-23				● Natural						
			75205-04		75206-04		75207-04				75208-04			● Black							
			73636-00		73636-01		73636-02				73636-03			● Diamond Pattern Black							
			75245-20		75246-20		75247-20				75248-20			● Natural							
			75245-04		75246-04		75247-04				75248-04			● Black							
	UltraSheer SENSITIVE Stockings Thigh High		77686-12		77686-13		77686-14				77686-15			15-20 mmHg	● Natural						
			77686-04		77686-05		77686-06				77686-07				● Black						
			77687-12		77687-13		77687-14				77687-15			20-30 mmHg	● Natural						
			77687-04		77687-05		77687-06				77687-07				● Black						
			77688-12		77688-13		77688-14				77688-15			30-40 mmHg	● Natural						
			77688-04		77688-05		77688-06				77688-07				● Black						
	UltraSheer Stockings Waist High		75173-12		75174-12		75175-12				75176-12			15-20 mmHg	● Natural						
			75173-04		75174-04		75175-04				75176-04				● Black						
		73637-00		73637-01		73637-02				73637-03			● Diamond Pattern Black								
		75209-13		75210-13		75211-13				75212-13			● Natural								
		75209-04		75210-04		75211-04				75212-04			● Black								
		73638-00		73638-01		73638-02				73638-03			● Diamond Pattern Black								
		75249-13		75250-13		75251-13				75252-13			● Natural								
		75249-04		75250-04		75251-04				75252-04			● Black								
UltraSheer Stockings Maternity		75177-04		75178-04		75179-04				75180-04			15-20 mmHg	● Natural							
		75177-01		75178-01		75179-01				75180-01				● Black							
		75213-06		75214-06		75215-06							20-30 mmHg	● Natural							
		75213-03		75214-03		75215-03								● Black							
DAILY	forMen Casual Socks Knee High		75476-02		75477-02		75478-02				75479-02			15-20 mmHg	● Black						
			75476-03		75477-03		75478-03				75479-03				● Khaki (Beige)						
			75484-02		75485-02		75486-02				75487-02			20-30 mmHg	● Black						
			75484-03		75485-03		75486-03				75487-03				● Khaki (Beige)						
			75492-02		75493-02		75494-02				75495-02			30-40 mmHg	● Black						
			75492-03		75493-03		75494-03				75495-03				● Khaki (Beige)						
WORK WEAR	forMen Socks Knee High		75253-01		75254-01		75255-01				75256-01			15-20 mmHg	● Black						
			75257-01		75258-01		75259-01				75260-01				● Black						
			75265-01		75266-01		75267-01				75268-01			30-40 mmHg	● Black						
			75261-00		75262-00		75263-00				75264-00				● Black						
forMen Socks Thigh High		75269-00		75270-00		75271-00				75272-00			30-40 mmHg	● Black							
		75147-03		75148-03		75149-03				75150-03			15-20 mmHg	● Black							
SPORT	ActiveWear Socks Knee High		75147-02		75148-02		75149-02				75150-02				20-30 mmHg	○ White					
			75155-03		75156-03		75157-03				75158-03			● Black							
			75155-02		75156-02		75157-02				75158-02			○ White							
			75411-03		75412-03		75413-03				75414-03			30-40 mmHg	● Black						
			75411-02		75412-02		75413-02				75414-02				○ White						
		TRADITIONAL	Relief Garments Knee High		75339-12		75340-12		75341-12		75343-12		75342-12		75344-12	15-20 mmHg	● Black				
	75339-10				75339-11		75340-10		75341-10		75343-10		75344-11	● Beige							
	78050-04				78050-15		78051-04		78051-15		78052-04		78052-15	77975-15	77975-16	78053-04	78053-15	77976-15	77976-16	20-30 mmHg	● Beige
	78054-04				78054-15		78055-04		78055-15		78056-04		78056-15	77977-15	77977-16	78057-04	78057-15	77978-15	77978-16		30-40 mmHg
Relief Garments Thigh High			78042-08		78042-09		78043-08		78043-09		78044-08		78044-09		78045-08	78045-09			20-30 mmHg	● Beige	
			78059-08		78059-09		78047-08		78047-09		78048-08		78048-09		78049-08	78049-09			30-40 mmHg	● Beige	
Relief Garments Waist High			78476-10		78476-11		78477-10		78477-11		78478-10		78478-11		78479-10	78479-11			20-30 mmHg	● Beige	
			78480-10		78481-10		78482-10				78483-10								30-40 mmHg	● Beige	
Relief Garments Chaps			78031-05		78032-05		78033-05				78034-05				20-30 mmHg	● Beige Left					
			78031-06		78032-06		78033-06				78034-06					● Beige Right					
			78035-05		78036-05		78037-05				78041-05				30-40 mmHg	● Beige Left					
			78035-06		78036-06		78037-06				78041-06					● Beige Right					

Selecting the product that serves your patient's needs




Different types of graduated compression socks or stockings are required for different indications.

The **JOBST® Medical LegWear range** is designed to relieve tired, aching legs, help treat venous conditions, and help prevent swelling of the feet, ankles and legs.

A measured amount of graduated compression is applied to your legs. It is highest at the ankle and tapers off up the leg, promoting normal flow of venous blood up the leg. It is therefore ideal for patients who are mobile and requiring support of the venous system to reduce oedema.

Alternatively, anti-embolic stockings are intended for use during surgical procedures and on patients confined to bed. They are designed to prevent blood clots rather than active oedema management.

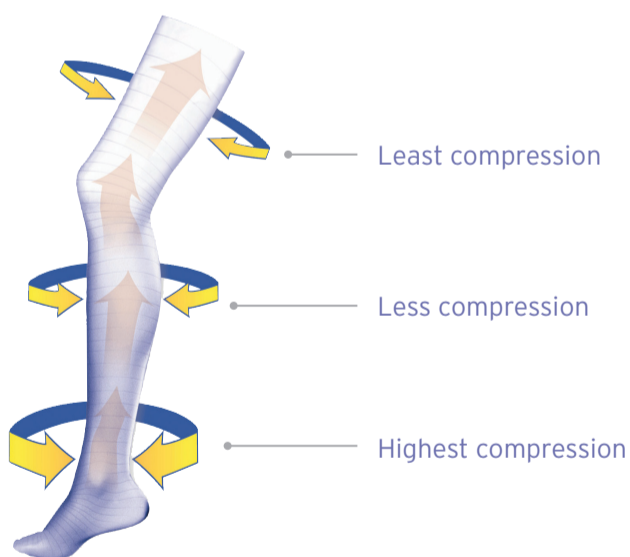

1 Select the required compression

15-20mmHg

- Tired aching legs
- Minor ankle, leg and foot swelling
- Minor varicose during pregnancy
- Long distance travel

20-30mmHg

- Tired aching legs
- Moderate to severe varicose during pregnancy
- Post vein surgery/post sclerotherapy
- Moderate tendency of oedema
- Superficial thrombophlebitis
- Helps prevent recurrence of ulcerations

30-40mmHg

- Severe varicose veins
- Lymphatic oedema
- Post vein surgery/post sclerotherapy
- Venous ulcers
- Manages manifestations of Post Thrombotic Syndrome (PTS), Chronic Venous Insufficiency CVI/CVD
- Aids in prevention of thrombotic syndrome

JOBST® Medical LegWear should not be worn and is contraindicated for the following conditions:

- Ischemia (eg. Advanced arterial disease)*
- Uncontrolled congestive heart failure
- Untreated septic phlebitis
- Phlegmasia coerulea dolens

*As part of a comprehensive assessment, a Doppler Ultrasound is used in the vascular screening process to clarify the vascular aetiology. A hand-held Doppler will assist in determining the degree, if any, of arterial disease. Using the Doppler the Ankle/Brachial Index (ABPI) can be calculated.

2 Select garment style



* open toe available


Some considerations when selecting the most appropriate style for your patient:

- Indication being treated/managed
- Lifestyle, mobility and age
- Ability to don and tolerate level of compression
- Skin Integrity
- Open toe or closed toe requirements
- Maternity Style only available in JOBST® UltraSheer

3 Measure for JOBST Medical LegWear

Measuring positions

- cH - Hip**
Circumference at largest dimension around buttocks.
- cG - Thigh**
Circumference at largest dimension.
- cC - Calf**
Circumference at largest dimension.
- cB - Ankle**
Circumference at smallest dimension.

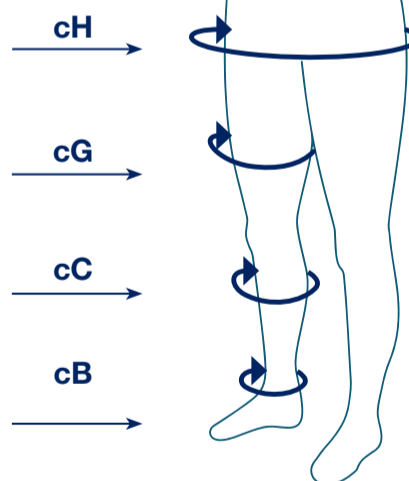

Provide the following measurements for:

- Knee High - cB & cC
- Thigh High - cB, cC & cG
- Waist high - cB, cG & cH

- Each leg is measured independently
- Measure in the early morning when possible as swelling occurs later in the day

JOBST UltraSheer, forMen, forMen Casual, ActiveWear and Relief Sizing

FC: Full Calf

Style	Measurement	Small	Medium	Large	Large FC	XLarge	XLarge FC
Unisex	Ankle	18 - 21cm	21 - 25cm	25 - 29cm	25 - 30cm	29 - 33cm	30 - 36cm
	Calf	28 - 38cm	30 - 42cm	32 - 46cm	46 - 61cm	34 - 50cm	46 - 61cm
	Thigh	40 - 62cm	46 - 70cm	54 - 78cm		60 - 81cm	
	Hip	71 - 117cm	76 - 127cm	81 - 137cm		102 - 166cm	
Maternity	Ankle	18 - 21cm	21 - 25cm	25 - 29cm		29 - 33cm	
	Thigh	40 - 62cm	46 - 70cm	54 - 78cm		60 - 81cm	
	Hip	81 - 132cm	86 - 142cm	91 - 152cm		112 - 180cm	